



KRUUSE Sacryl Plus

Multifilament, braided: For secure knot tying and easy handling

Synthetic: Minimal tissue reaction

Coated: For easy knot tying and smooth passage through tissue

Antibacterial: To reduce incidence of infection, oedema and pain

Absorbable: By hydrolysis

KRUUSE Sacryl Plus is an antibacterial, absorbable, synthetic, braided and coated suture containing Chlorhexidine diacetate. Chlorhexidine diacetate is the antibacterial component.

KRUUSE Sacryl Plus, being an absorbable suture, should not be used where extended approximation of tissue under stress or tension is required and in patients with known allergic reactions to Chlorhexidine diacetate.

KRUUSE Sacryl Plus tested in vitro creates a zone of inhibition against the most common surgical site pathogens including Staphylococcus aureus.

Absorption: Absorbable:

- Complete mass absorption in approximately 56-63 days
- Absorption by hydrolysis

Chemical composition:

- KRUUSE Sacryl is composed of Polyglactin 910 (copolymer of 90% Glycolide and 10% L-lactide)

Tensile Strength:

- 70 % remains after 14 days
- 50 % remains after 21 days
- 25 % remains after 28 days
- Complete mass absorption in approximately 56-63 days

Colour:

- Violet and undyed

Indications:

KRUUSE Sacryl sutures are indicated for soft tissue approximation and ligation, where absorbable sutures are indicated, including:

- Suturing of subcutaneous tissue
- Fascia
- Joint capsules
- Muscles
- Uterus and GI-surgery
- Vascular surgery

Contraindications:

- Where permanent support is required

Range:

- USP 5-0 to 2
- EP 1 to 5

Sterilization method:

- EO

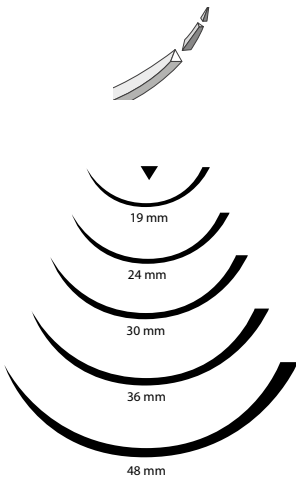
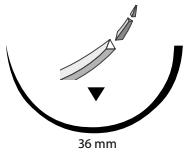

Packing:

- Boxes of 18 peel packs

v



KRUUSE Sacryl Plus

Shape of needle	Needle length mm	USP	EP	Length in cm	Suture colour	Number per box	Cat. No
<p>3/8 circle, reverse cutting</p> 	19	5-0	1	45	Undyed	18	152339
	19	4-0	1.5	45	Undyed	18	152338
	19	4-0	1.5	70	Violet	18	152336
	19	3-0	2	70	Violet	18	152333
	19	2-0	3	70	Violet	18	152340
	24	3-0	2	70	Undyed	18	152337
	24	3-0	2	70	Violet	18	152334
	24	2-0	3	70	Violet	18	152331
	30	3-0	2	70	Violet	18	152335
	30	2-0	3	70	Violet	18	152332
	30	0	3.5	70	Violet	18	152330
	36	1	4	90	Violet	18	152341
	48	2	5	90	Violet	18	152342
	<p>1/2 circle, reverse cutting</p> 	36	0	3.5	70	Violet	18
36		1	4	70	Violet	18	152346
<p>1/2 circle, round bodied, tapercut</p> 	26	3-0	2	70	Violet	18	152343
	26	2-0	3	70	Violet	18	152344