

BUSTER Foley Plus Catheter

2-way catheter with retention balloon, high biocompatibility and suitable for long term placement.

Perfect for

- Drainage purposes
- Urine sampling
- Bladder lavage

User-friendly

- Easy placement
- Radiopaque
- Size colour-coding

Patient comfort

- Soft and flexible
- High biocompatibility
- Long-term placement


Improved patient safety

The retention balloon secures the catheter in situ and the rounded tip and atraumatic side holes help avoid tissue damage. The highly biocompatible silicone makes the catheter suitable for long-term placement up to 7 days.


Easy insertion - no wire

The catheter is made of white silicone with an added substance that increases the shore of the catheter. This means that the catheter rigidity is increased and therefore does not require a wire for insertion.




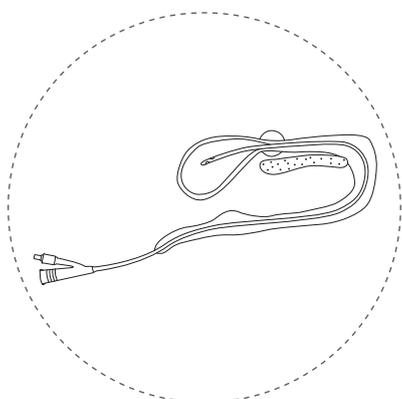

Insert the catheter through urethra into the bladder

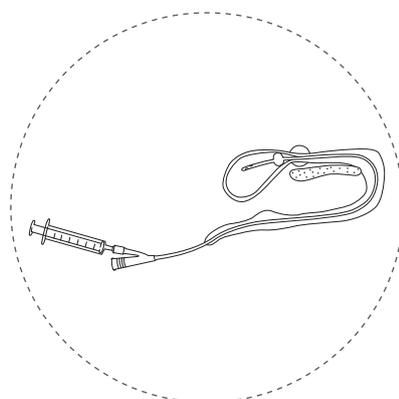

Fill the balloon with isotonic saline

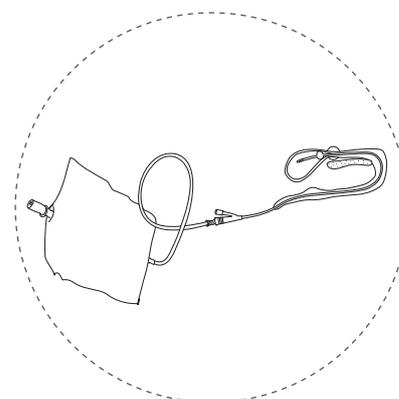

Connect a BUSTER Urine Bag to collect and control urine flow

Cat. No	Description	Type	Size, Metric	Size, Imperial	Colour	Package
273880	BUSTER Foley Plus Catheter	Side holes, without wire	2 mm x 30 cm	6 Fr x 12 in	●	5/pk
273881	BUSTER Foley Plus Catheter	Side holes, without wire	2 mm x 55 cm	6 Fr x 22 in	●	5/pk
273882	BUSTER Foley Plus Catheter	Side holes, without wire	2.7 mm x 30 cm	8 Fr x 12 in	●	5/pk
273883	BUSTER Foley Plus Catheter	Side holes, without wire	2.7 mm x 55 cm	8 Fr x 22 in	●	5/pk
273884	BUSTER Foley Plus Catheter	Side holes, without wire	3.3 mm x 30 cm	10 Fr x 12 in	●	5/pk
273885	BUSTER Foley Plus Catheter	Side holes, without wire	3.3 mm x 55 cm	10 Fr x 22 in	●	5/pk

