

## BUSTER Dog Flushing Catheter, Open End

Based on the popular BUSTER Dog Catheter, the new open-end catheter is ideal for flushing of the urethra to relieve urethral obstruction.

### Perfect for

- Flushing of the urethra
- Passage check
- Bladder lavage
- Urine sampling

### Helps relieve

- Urethral obstruction
- Urinary retention

### Minimal patient discomfort

- Easy placement in situ
- Minimal pressure on the urethral wall during flushing

### Best practice

An open-end catheter provides optimal conditions for flushing uroliths retrograde into the bladder.

### User-friendly

The polypropylene catheter tubing is rigid and still flexible for easy catheter introduction and placement. The female Luer Lock fitting ensures easy and secure syringe attachment.

### Increased patient safety

This catheter helps prevent unnecessary tissue damage thanks to the atraumatic, rounded end hole and because the hydro-pressure is directly aimed at the obstruction, not the urethral wall.




Andrew Jackson, DVM, DACVS states:

“Urethral obstructions in male dogs are common in veterinary practice. These obstructions can be difficult to relieve, and many patients are referred to emergency/specialty hospitals still obstructed. It doesn’t have to be this way! With a few precise but straightforward steps, general practitioners can relieve most urethral obstructions with minimal fuss to the patient.”


Source: <https://www.dvm360.com/view/surgery-stat-pinch-and-flush-relieving-urethral-obstructions-male-dogs>




Atraumatic and polished open-end hole for optimal and correct flushing direction

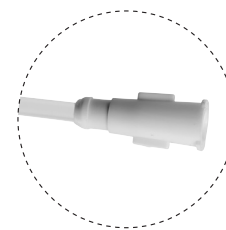

Female Luer Lock fitting for easy and secure syringe attachment

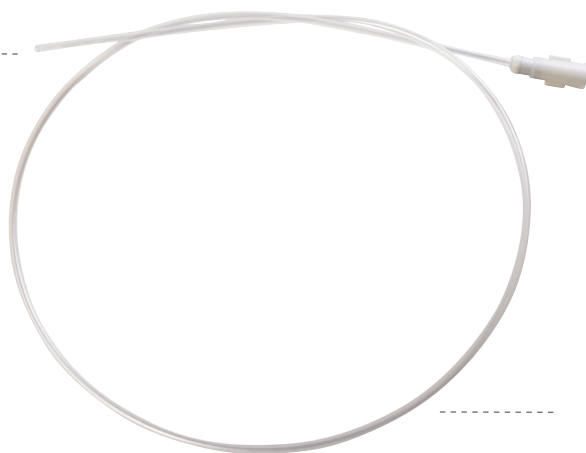

Rigid and flexible PP tubing for easy introduction and placement

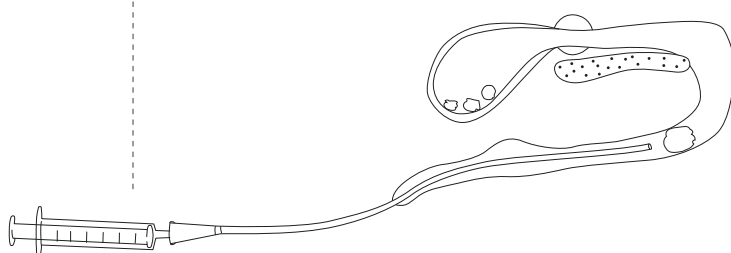

Cat. No	Description	Type	Size, Metric	Size, Imperial	Package
273415	BUSTER Dog Flushing Catheter	Open End	1.3 mm x 50 cm	4 Fr x 20 in	12/pk
273416	BUSTER Dog Flushing Catheter	Open End	2.0 mm x 50 cm	6 Fr x 20 in	12/pk
273417	BUSTER Dog Flushing Catheter	Open End	2.6 mm x 50 cm	8 Fr x 20 in	12/pk
273418	BUSTER Dog Flushing Catheter	Open End	3.3 mm x 50 cm	10 Fr x 20 in	12/pk