



## KRUUSE Nylon

**Monofilament:** For smooth passage through tissue and high tissue compatibility

**Synthetic:** Minimal tissue reaction

**Non-absorbable:** Provides permanent tensile strength

KRUUSE Nylon sutures are coloured blue for easy identification during surgery. The monofilament and synthetic properties of the suture give a smooth surface, minimal tissue reaction and decreased incidence of infection and trauma.

KRUUSE Nylon suture is quite pliable and does not support bacterial growth due to the smooth surface. The elasticity of KRUUSE Nylon suture makes it useful in retention and skin closure.

**Absorption:**

- Non-absorbable

**Chemical composition:**

- KRUUSE Nylon monofilament suture is made of monofilament polyamide 6&66 by extruding a strong uniform strand.

**Colour:**

- Blue

**Indications:**

KRUUSE Nylon is commonly used in soft tissue approximation and ligation where permanent support is needed. Areas of use include:

- Skin closure
- Ophthalmic surgery
- Vascular and cardiovascular surgery
- Neurosurgery

**Range:**

- USP 5-0 to 2
- EP 1 to 5

**Sterilization method:**

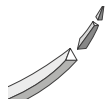
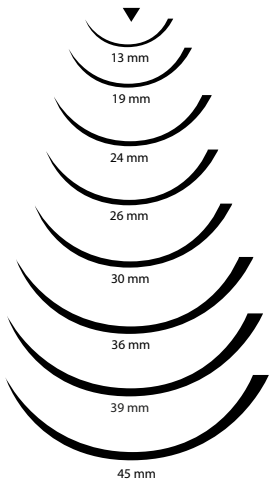

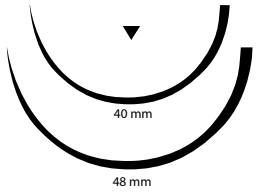


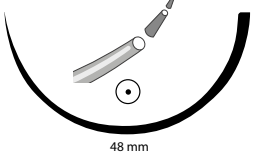
- EO

**Packing:**

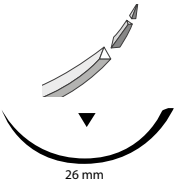
- Boxes of 18 peel packs



## KRUUSE Nylon

Shape of needle	Needle length mm	USP	EP	Length in cm	Suture colour	Number per box	Cat. No
<p>3/8 circle, reverse cutting</p>  	13	5-0	1	45	Blue	18	152033
	19	5-0	1	45	Blue	18	152417
	19	4-0	1.5	45	Blue	18	152030
	19	3-0	2	45	Blue	18	152031
	19	3-0	2	70	Blue	18	152507
	24	4-0	1.5	70	Blue	18	152425
	24	3-0	2	70	Blue	18	152508
	24	0	3.5	70	Blue	18	152420
	26	3-0	2	70	Blue	18	152424
	26	2-0	3	70	Blue	18	152423
	30	2-0	3	90	Blue	18	152429
	30	0	3.5	70	Blue	18	152421
	36	1	4	70	Blue	18	152422
	39	2-0	3	70	Blue	18	152418
45	2-0	3	90	Blue	18	152032	
<p>1/2 circle, reverse cutting</p>  	40	0	3.5	90	Blue	18	152426
	48	1	4	90	Blue	18	152427
<p>3/8 circle, conventional cutting</p> 	90	1	4	100	Blue	18	152730
	90	2	5	90	Blue	18	152428
	90	2	5	200	Blue	18	152731
<p>Straight, cutting</p> 	55	2-0	3	90	Blue	18	152419
<p>1/2 circle, round bodied, taperpoint</p> 	48	2	5	90	Blue	18	152509

## KRUUSE Nylon - Race Track

Shape of needle	Needle length mm	USP	EP	Length in cm	Suture colour	Number per box	Cat. No
3/8 circle, reverse cutting 	26	3-0	2	70	Blue	15	152840
	26	2-0	3	70	Blue	15	152841

